


Basis Q

trifocal

Single piece, sterile, foldable acrylic aspheric hydrophilic posterior chamber intraocular lens (IOL) with UV blocker and blue light filter. The optic is trifocal, with 3 mm diameter diffractive array on the anterior surface utilizing EPS technology. The IOL is non-preloaded, which has to be manually loaded into a compatible injector before implantation.

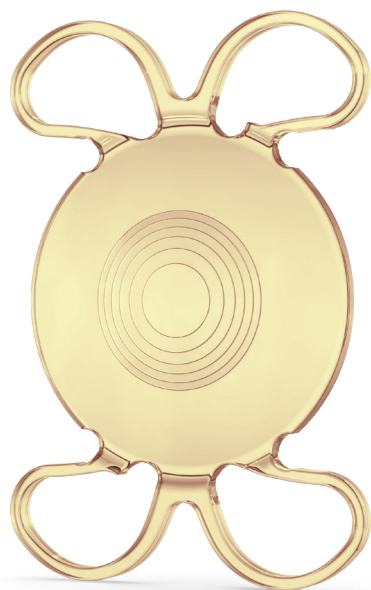

MATERIAL

Copolymer of hydrophobic and hydrophilic monomers with 25% water content, UV blocker and blue light filter

Refractive Index

1.46

Abbe Number

58

DESIGN

Optic Diameter
Overall Diameter

6 mm

11.0 mm (8.0 D → +15.0 D)
10.7 mm (+15.5 D → +35.0 D)

Optic Design

Aspheric - Aberration Neutral
Biconvex (+8.0 D → +35.0 D)

Nominal A-constant

118,9

Haptic Design

Posterior vaulting 4 closed loops
with 0° angulation

PCO Prevention

360° sharp edge

Shelf Life

5 years after sterilization

INJECTION

Estimated Incision Size

1.8 - 2.2 mm

Recommended Injector

Compatibility chart available:
1stq.eu/compatibility/

CONSTANTS

for optical and immersion ultrasound measurements

SRK/T A-Constant

A = 118.9

Haigis

a0 = 1.243
a1 = 0.400
a2 = 0.100

Hoffer Q

pACD = 5.46

Holladay I

SF = 1.67

Barrett Universal II


LF = 1.83

Holladay II

ACD = 5.49

Basis Q

trifocal

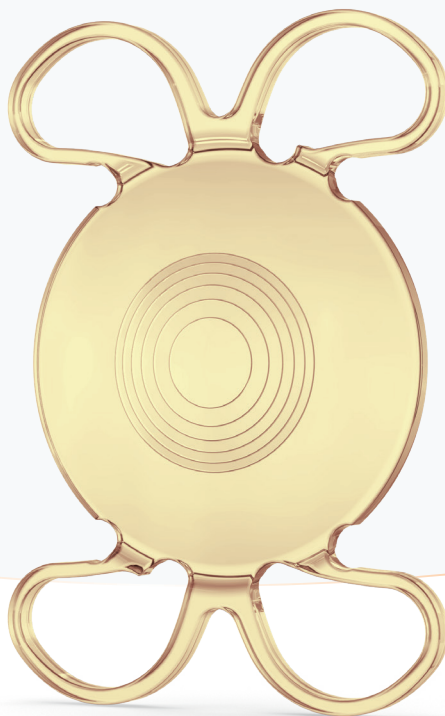

POWER RANGE

Article No. Yellow	SEQ/Powers	Steps	Cylinders	Steps	Addition
S2EWYN	+8.0 D to +9.0 D	1.0 D			+3.5 D near +1.75 D intermediate
B2EWYN	+10.0 D to +30.0 D	0.5 D			+3.5 D near +1.75 D intermediate
S2EWYN	+31.0 D to +35.0 D	1.0 D			+3.5 D near +1.75 D intermediate



1stQ GmbH

Konrad-Zuse-Ring 23 | 68163 Mannheim
Germany



Leading
IOL solutions

www.1stq.eu
www.1stq.de/en